

Product datasheet

Rabbit Anti Streptavidin Antibody, pAb

Cat: NA0077

Name: Rabbit Anti Streptavidin Antibody, pAb

Immunogen: Recombinant streptavidin

Host Species: Rabbit

Specificity: This antibody is specific to streptavidin.

Form: Liquid

Concentration: 1 mg/ml

Storage Buffer: 1×PBS, pH7.4, containing 50% glycerol and 0.03% Proclin 300.

Purification: Antigen affinity purification

Clonality: Polyclonal

Isotype: Rabbit IgG

Shipping and Storage: Shipped in ice bag. The antibody can be stored for 1-2 weeks at 2-8°C. For long term storage, store it at -20°C. Avoid repeated freezing and thawing cycles.

Expiration Date: Specified on product label.

Applications

Application	Recommended Usage
ELISA	0.05-0.5µg/ml
Western Blot	0.1-1µg/ml
Other applications	User-optimized

Note: Ningyang Biotech can offer this product according to your requirement, including product size, buffer components, etc.

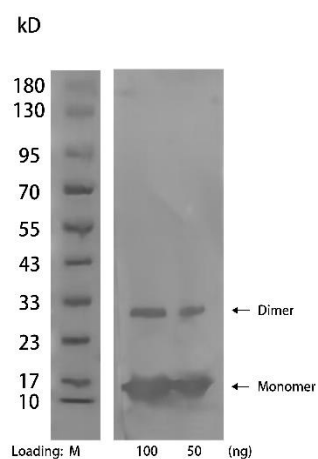


Figure 1. Detection of streptavidin fusion protein were analyzed by Western blot using Rabbit Anti Streptavidin Antibody, pAb (NingyangBiotech,Cat.NA0077, at 1 μ g/mL concentration).

The signal was developed with Goat Anti-Rabbit IgG(H&L)[HRP], pAb(NingyangBiotech Cat.NA0030).

For research use only. Not intended for human and animal therapeutic or diagnostic use.