

Product datasheet

Rabbit Anti T785(AXC-715) Antibody, mAb

Cat: NA0178

Name: Rabbit Anti T785(AXC-715) Antibody, mAb

Immunogen: T785(AXC-715) conjugated to keyhole limpet haemocyanin(KLH)

Host Species: Rabbit

Specificity: This antibody recognizes T785(AXC-715) and T785(AXC-715) conjugate.

Form: Liquid

Concentration: 0.5 mg/ml

Storage Buffer: 1×PBS, pH7.4, containing 0.1% BSA, 50% glycerol and 0.03% Proclin 300.

Purification: Purified by Protein A

Clonality: Monoclonal

Clone ID: 11C6

Shipping and Storage: Shipped in ice bag. The antibody can be stored for 1-2 weeks at 2-8°C. For long term storage, store it at -20°C. Avoid repeated freezing and thawing cycles.


Expiration Date: Specified on product label.

Applications

Application	Recommended Usage
ELISA	0.01-0.1 µg/ml
Other applications	User-optimized

Note: Ningyang Biotech can offer this product according to your requirement, including product size, buffer components, etc.

Examples

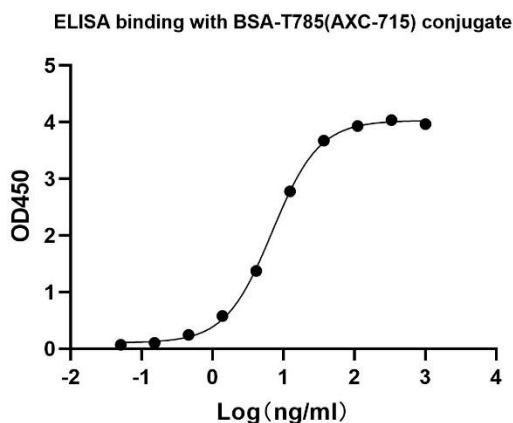


Rabbit Anti T785(AXC-715) Antibody, mAb (Cat:NA0178) binds with BSA-T785(AXC-715) conjugate.

Coating antigen: BSA-T785(AXC-715) conjugate, 1 µg/ml, 100 µL/well.

Rabbit Anti T785(AXC-715) Antibody, mAb (Cat:NA0178) dilution start from 1000 ng/ml.

EC₅₀= 6.92 ng/ml.

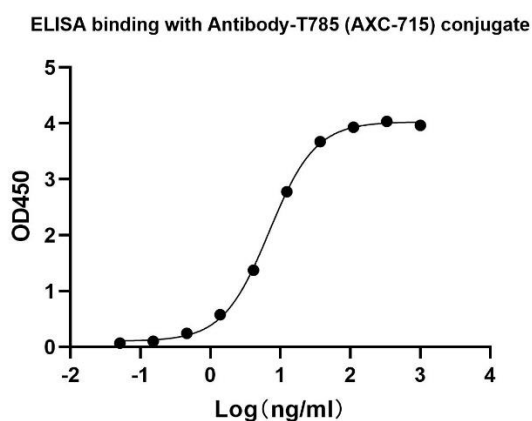

Rabbit Anti T785(AXC-715) Antibody, mAb (Cat:NA0178) binds with Antibody-T785(AXC-715) conjugate.

Coating antigen: Antibody-T785(AXC-715) conjugate, 1 µg/ml, 100 µL/well.

Rabbit Anti T785(AXC-715) Antibody, mAb (Cat:NA0178) dilution start from 1000 ng/ml.

EC₅₀= 2.97 ng/ml.

For research use only. Not intended for human and animal therapeutic or diagnostic use.